


carilex®

CariChair

A semi-automatic mobile cushion system with battery backup defies time and place constraints for uninterrupted, around-the-clock care. Whenever, wherever. With CariChair, the patients no longer need to be confined to the bed for pressure ulcer management. An exceptionally progressive accessory to the extensive Carilex wound care solution family, CariChair is designed to be compatible with various mobile devices, including geriatric chairs, recliners, and wheelchairs.

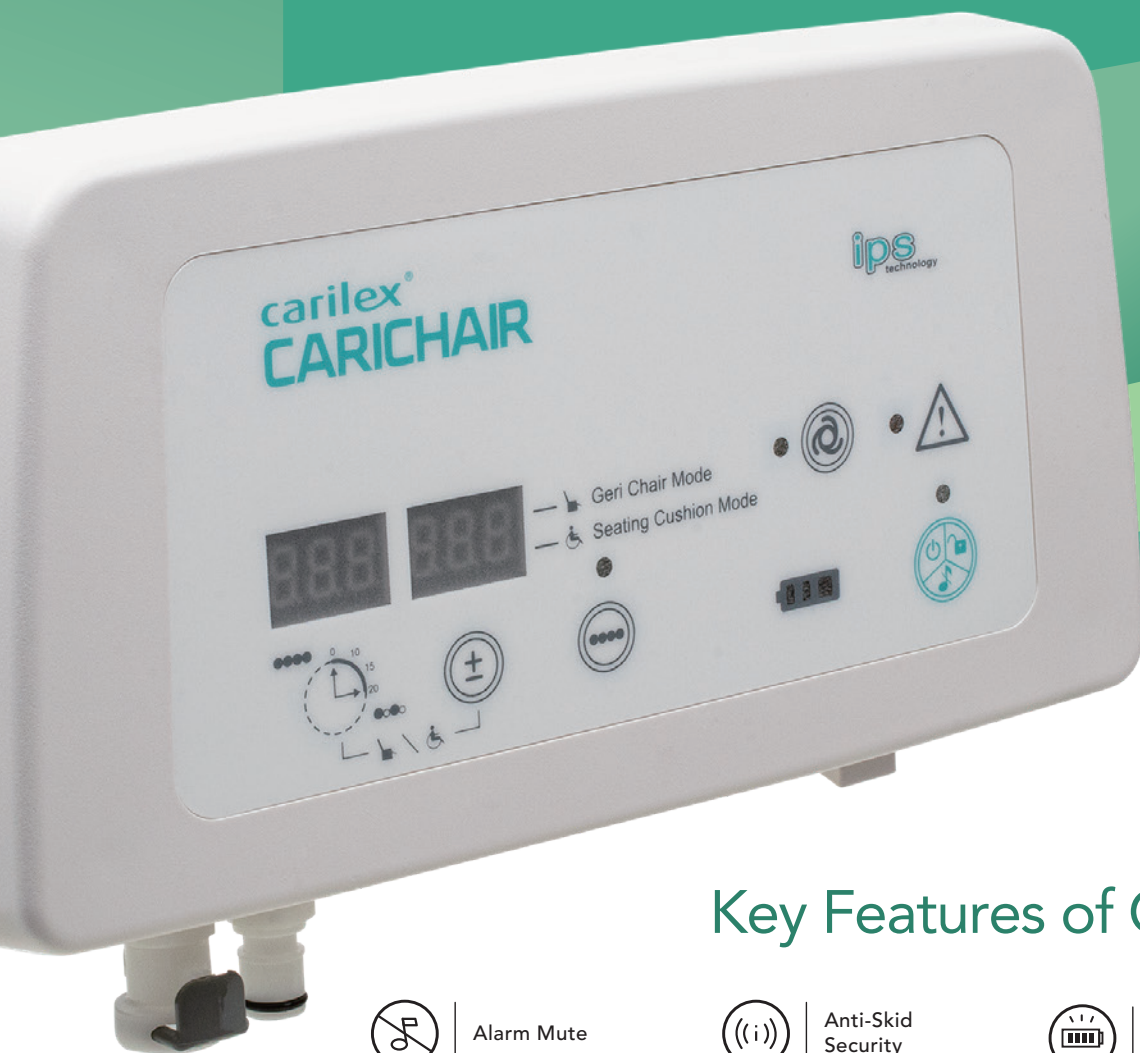

Key Features of CariChair



Alarm Mute



Anti-Skid Security



Battery Level LED



Carry bag



Digital LED Screen



Panel Lock



Visual and Audible Indicator



Whisper Quiet



Washable, Waterproof Cushion Cover

Multiple Compatibility

Compatible with multiple cushions for various types of chairs

Intelligent Pressure Sensing

Responds to patient movements on cushions by automatically adjusting internal mattress pressure

MobiCare Plus C-set



Battery Backup Up to 8 Hours

Automatic power detection sensor activates battery backup, in the event of main power failure

Dynamic Mode

Alternating therapy

Static Therapy/CLP

All air cells maintain constant pressure support

Comfort Setting

Provides optimal internal Pressure according to each patient's sensation. Includes options for comfort level adjustment, no weight input needed

Multiple Cycle Times

10, 15 and 20 minute cycles for alternating pressure therapy

Technical Data	Power Unit Model	REF : SR321
	Power Input	100-240Vac 50-60Hz 1.5A
	Power Consumption	23W Adaptor + Pump; 7W Pump only
	Weight	Power Unit : 0.92Kg (2.0lb) Cushion Set : • CariChair Cushion 2.1Kg (4.7lb) • MobiCare Plus / MobiCare Cushion 1.5Kg (3.3lb)
	Dimension	Power Unit : 10" x 6" x 2" (24.5 x 15.5 x 6.1cm) (W x H x D) Cushion Set : (W x L x H) • CariChair Cushion, 71" x 20" x 6" (180 x 50 x 14 cm) • MobiCare Plus Cushion, 18" x 18" x 4" (45 x 45 x 9 cm) • MobiCare Cushion, 18" x 18" x 4" (45 x 45 x 9 cm)
	Agency Approval	US NRTL - MET classified Mark_ANSI/AAMI ES60601-1
	Operating Condition	• Temperature : 5°C to 40°C (41°F to 104°F) • Humidity : 15% to 90%
	Transportation & Storage Condition	• Temperature : -25°C to 70°C (-13°F to 158°F) • Humidity : 0% to 90%
	Weight Capacity	• CariChair PU+CariChair C-set: 150Kg (330lb) • CariChair PU+MobiChair Plus set: 150Kg (330lb) • CariChair PU+MobiCare C-set: 120Kg (260lb)

This Power Unit is not AP/APG protected