

# VT Dressing Kits

Carilex Dressing Kits are indicated for the patient who benefits from wound management therapy via the application of negative pressure wound therapy to remove fluids and excess exudates, infectious material, and tissue debris which may promote wound healing.

The VT Dressing Kits have been carefully engineered and tested for comfort and safety while delivering superior performance by maintaining prescribed negative pressure at the wound area.



## VT Port Pad

- ▶ The VT port pad is made of soft medical-grade PVC
- ▶ The port comes with a circle drape and customized port seal for easy application
- ▶ The design of the port minimizes shear forces on the peri-wound and maximizes flow of exudate from the wound while maintaining a low profile for greater patient comfort
- ▶ The pad is flexible and readily contours to the application site on the body



## VT Foam

- ▶ The VT Foam is reticulated hydrophobic, polyurethane foam that aids the formation of granulation tissue due to its cell structure and pore size



## VT Drape

- ▶ The VT drape is manufactured from a high quality film and is specifically designed for application in NPWT

**Carilex VT Dressing Kits are available with options for 3 different sizes**

### Large



### Medium



### Small



## ORDERING INFORMATION

PRODUCT NUMBER	PRODUCT NAME	DESCRIPTION
S1001-2040	VT Dressing Kit S Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size small(4" x 3.2" x 1.2" (100 x 80 x 30mm)) with lockable connector.</li> <li>1 port set and 1 pc. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2050	VT Dressing Kit M Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size medium(7.9" x 5" x 1.2" (20 x 125 x 30mm)) with lockable connector.</li> <li>1 port set and 2 pcs. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2060	VT Dressing Kit L Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size large (9.9" x 5.9" x 1.2" (250 x 150 x 30mm)) with lockable connector.</li> <li>1 port set and 3 pcs. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2100	VT Spiral Dressing Kit Black Lock	<ul style="list-style-type: none"> <li>2 pcs. of black foam dressing kits with spiral(5.7" x 0.6" (146 x 15mm)) patent cutting &amp; lockable connector.</li> <li>1 port set and 2 pcs. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2120	VT Dressing Kit XL Carilex	<ul style="list-style-type: none"> <li>2 pcs. of black foam dressing kits, size extra large(9.6" x 14.4" x 0.6" (245 x 365 x 15mm)) with lockable connector.</li> <li>1 port set and 3 pcs. of 11.9" x 15.8" (300 x 400 mm) drape</li> </ul>
S1001-2130	VT Drape and Port Set Kit S Lock	<ul style="list-style-type: none"> <li>1 port set and 1 pc. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2140	VT Drape and Port Set Kit M Lock	<ul style="list-style-type: none"> <li>1 port set and 2 pcs. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2150	VT Drape and Port Set Kit L Lock	<ul style="list-style-type: none"> <li>1 port set and 3 pcs. of 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>
S1001-2200	VT Dressing Kit TS S Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size small with lockable connector.</li> <li>1 port set and 1 pc. of 7.9" x 11.9" (200 x 300 mm) drape, 2 step circle drape</li> </ul>
S1001-2210	VT Dressing Kit TS M Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size medium with lockable connector.</li> <li>1 port set and 2 pcs. of 7.9" x 11.9" (200 x 300 mm) drape, 2 step circle drape</li> </ul>
S1001-2220	VT Dressing Kit TS L Carilex	<ul style="list-style-type: none"> <li>Black foam dressing kits, size large with lockable connector.</li> <li>1 port set and 3 pcs. of 7.9" x 11.9" (200 x 300 mm) drape, 2 step circle drape</li> </ul>
S1001-8160	VT Port Set Kit	<ul style="list-style-type: none"> <li>one port set kit</li> </ul>
S1001-8190	VT Instill Port Set Kit	<ul style="list-style-type: none"> <li>one Instill port set kit</li> </ul>
S1001-8200	VT · 100 Y - Connector	<ul style="list-style-type: none"> <li>A "Y" connector to use for all VT series to be able to work with two wounds</li> </ul>
S1001-8300	VT Drape 2pcs pack	<ul style="list-style-type: none"> <li>A sterilized pack of 2 pcs. 7.9" x 11.9" (200 x 300 mm) drape</li> </ul>

### ① VT Foam

Polyurethane (PU)

### ② VT Drape & Circle Drape

Polyurethane Film Coated with Acrylic Adhesive

### ③ Tube

NON-P PVC

### ④ VT Port


NON-P PVC

### ⑤ Clip

Polypropylene Hom(PP)

### ⑥ Connector

ABS

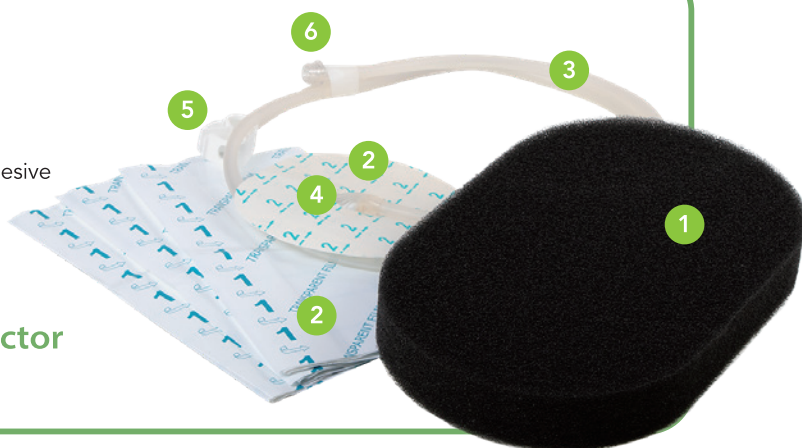

[www.carilexmedical.com](http://www.carilexmedical.com)

Carilex Medical, Inc. **HQ**

TEL : +886-3-328-7882  
Email : [ameo@carilexmedical.com](mailto:ameo@carilexmedical.com)  
No.77, Keji 1<sup>st</sup> Rd., Guishan District,  
Taoyuan City (33383), Taiwan

Carilex Medical USA

TEL : +1-866-551-4227  
Email : [americas@carilexmedicalusa.com](mailto:americas@carilexmedicalusa.com)  
17068 Evergreen Place, Unit A, City of Industry,  
CA 91745 USA

Carilex Medical DE

TEL : +49-4121-476896  
Email : [eu@carilexmedical.com](mailto:eu@carilexmedical.com)  
Fritz-Strassmann Str. 5, 25337 Elmshorn  
Germany

Carilex Medical NL

TEL : +31-2089-05788  
Email : [nl@carilexmedical.com](mailto:nl@carilexmedical.com)  
Overschiestraat 59, 1062XD  
Amsterdam

Carilex Medical UK

TEL : +44-1323-761900  
Email : [uk@carilexmedical.com](mailto:uk@carilexmedical.com)  
5<sup>th</sup> Floor, 167-169 Great Portland St,  
London W1W 5PF, United Kingdom